


# Подготовка к семинару по использованию Ion Reporter

## Создание аккаунта

1. На странице <https://ionreporter.thermofisher.com/ir/> нажмите “Register new account”

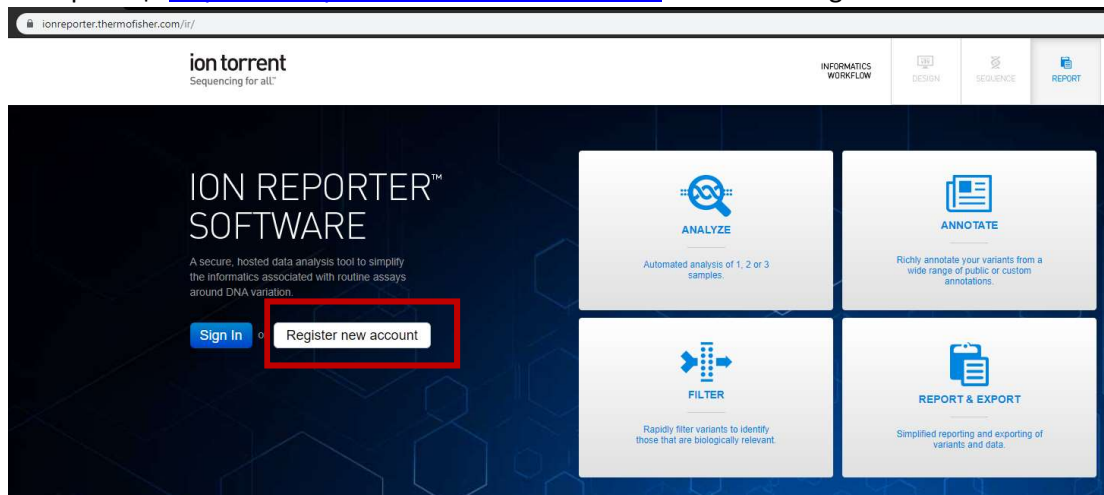

2. Заполните анкету и нажмите “Create account”

**ThermoFisher**  
SCIENTIFIC

Why you should create an account

Be the first to hear about events and training in your region, get access to the latest updates and product information, get full access to the Ion Reporter™, Ion Ampliseq™ Designer and Cloud Suite applications.

ThermoFisher SCIENTIFIC applied biosystems invitrogen gibco iontorrent

After creating a basic account, be sure to create a full account to get all these great features.

First name \*

Last name \*

Email (This will be your username) \*

Confirm your email \*

Password \*

Reenter password \*

Password strength:

Company or institution \*

Department / lab \*

City \*

Country/Region \*

Security question, if you forget your password \*

Security question answer \*

3. Вернитесь на <https://ionreporter.thermofisher.com/ir/> и создайте организацию  
Согласитесь с лицензионным соглашением:

Thermo Fisher Cloud Terms of Use

LIFE TECHNOLOGIES CORPORATION  
THERMO FISHER CLOUD  
TERMS OF USE  
Last Updated December 12, 2018

These Terms of Use ("Terms") govern your access to and use of Life's proprietary infrastructure, software, applications and services ("Services"). By accessing or using the Services, you accept these Terms and conclude a legally binding contract between you and Life Technologies Corporation, a Thermo Fisher Scientific business ("Life"). To access or use the Services, you must be 18 years or older and have the requisite power and authority to enter into these Terms, including on behalf of your organization if you use the Services for its benefit. If you do not accept these Terms, do not register for or use the Services.

If you are located within the People's Republic of China (including Hong Kong), a separate set of Terms of Use apply to you. Any claim, dispute or controversy of whatever nature between you and Life arising out of or relating to these Terms or the Services will be resolved by final and binding arbitration in accordance with Section 16 below. Please print these Terms for your records.

1. DEFINITIONS

A. Parties

"You" and "your" refer to you, the user of the Services. A "user" is someone who registers for, accesses, browses, crawls, scrapes, or in any way uses the Services. If your organization has authorized or otherwise permits you to access or use the Services for its benefit, "you" also includes your organization.

"We", "us" and "our" refer to Life.

B. Services

[Printer-friendly version](#) [Accept](#)

Согласитесь с созданием организации:

Confirm Organization creation

Ion Reporter requires each user to be a member of a new or an existing Organization. If you wish to become a member of an existing Organization, please contact an administrator of that existing Organization. If you wish to create a new Organization and be its administrator, please create an Organization on Thermo Fisher cloud. Click on 'Yes' to be redirected to the Organization creation page. Click 'No' to login as a different user.

[Yes](#) [No](#)

Введите название организации:

app.thermofisher.com/apps/um/#neworg

User Manager

Create a New Organization

Enter a name for your new organization. Your organization name must be unique, and should describe your organization to new team members.

[Cancel](#) [Create](#)


Уже можно пользоваться аккаунтом и зайти в Ion Reporter. Если захотите пригласить коллег в свою организацию, можно это сделать через action -> invite user

app.thermofisher.com/apps/um/ir/

User Manager

Email Address	First Name	Last Name	Status	Last Login	Action
Ivan.Ivanov2@gmail.com	Ivan	Ivanov	Enabled	2019-10-15T11:00:48.756Z	<a href="#">Invite User</a> <a href="#">User Details</a>

4. Зайдите в Ion Reporter <https://ionreporter.thermofisher.com/ir/>

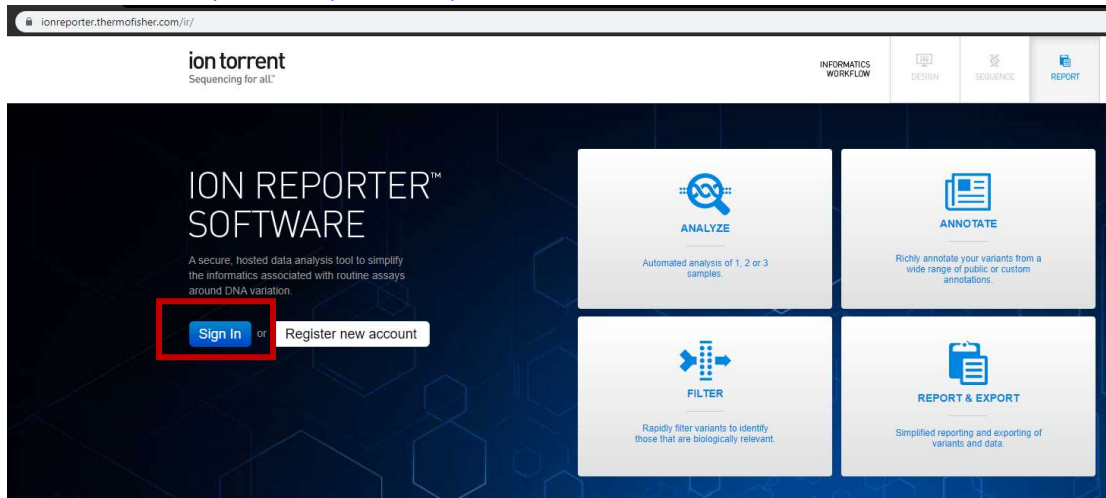

Согласитесь с лицензионным соглашением:

**License Agreement**

LIFE TECHNOLOGIES CORPORATION  
ION REPORTER™ SOFTWARE, VERSION 5.12.0  
TERMS OF USE  
Effective August 8, 2019

1. GENERAL

These terms of use ("Terms") are entered into by and between Life Technologies Corporation, a Thermo Fisher Scientific business ("Life", "we", "us" or "our") and you, the person registering for access to our proprietary infrastructure, software, applications and services (the "Platform"), or the organization or entity that has authorized you to register for access to the Platform for its benefit (in either case "you", "your"), and govern your access to and use of the Platform as a web-based hosted service made available by us (the "Service").


PLEASE READ THESE TERMS CAREFULLY. IT IS IMPORTANT TO UNDERSTAND THAT THIS IS A LEGALLY BINDING AGREEMENT BETWEEN YOU AND LIFE. BY REGISTERING FOR, ACCESSING or using THE SERVICE, YOU CONFIRM THAT (1) YOU HAVE READ THESE TERMS, (2) YOU ARE AT LEAST 18 YEARS OLD, (3) YOU, AND THE ORGANIZATION OR ENTITY ON WHOSE BEHALF YOU ARE ACCEPTING THESE TERMS AND/OR USING THE SERVICE, ACCEPT THESE TERMS, and (4) WHEN REGISTERING FOR THE SERVICE, YOU READ, UNDERSTOOD AND ACCEPTED OUR PRIVACY POLICY (available at <http://www.thermofisher.com/us/en/home/global/privacy-policy.html?cid=fl-privacy-policy>).

("Privacy Policy"). CLICKING ON THE "I ACCEPT" BUTTON BELOW INDICATES THAT YOU CONFIRM EACH OF 1-4 ABOVE, AND ACCEPT AND INTEND TO BE BOUND BY THESE TERMS, AND CONSENT TO ENTER INTO THESE TERMS IN ELECTRONIC FORM. IF YOU DO NOT


После принятия лицензионного соглашения вы окажетесь на домашней странице Ion Reporter.

## Запуск анализа


1. На домашней странице нажмите на вкладку “Analyses”

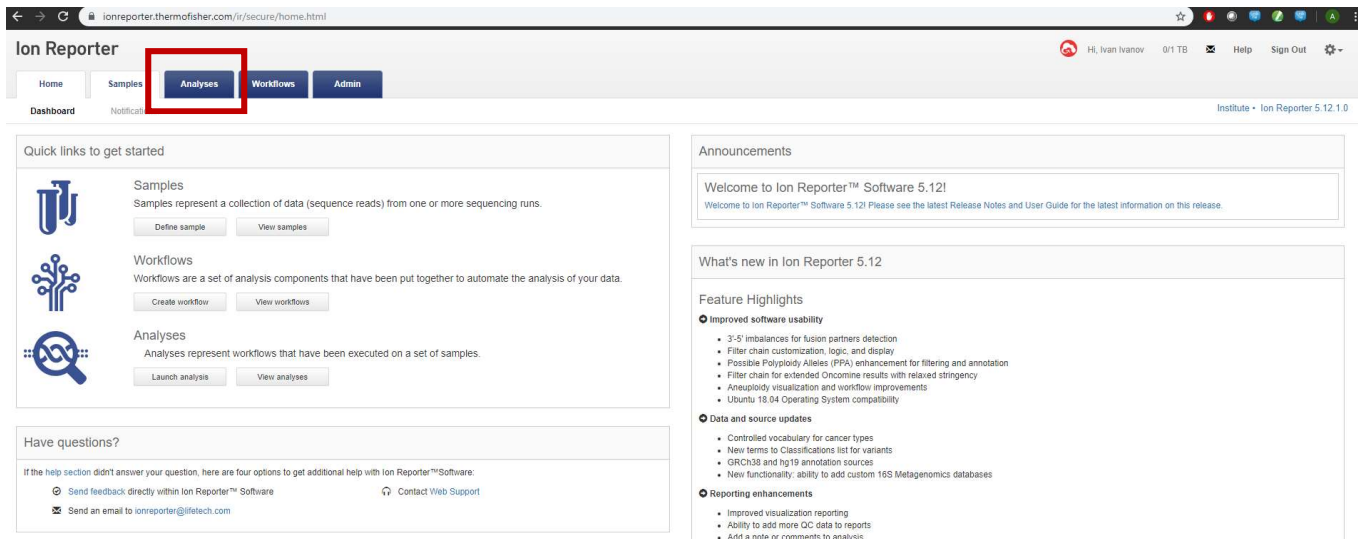

2. Во втором уровне вкладок нажмите “Launch”

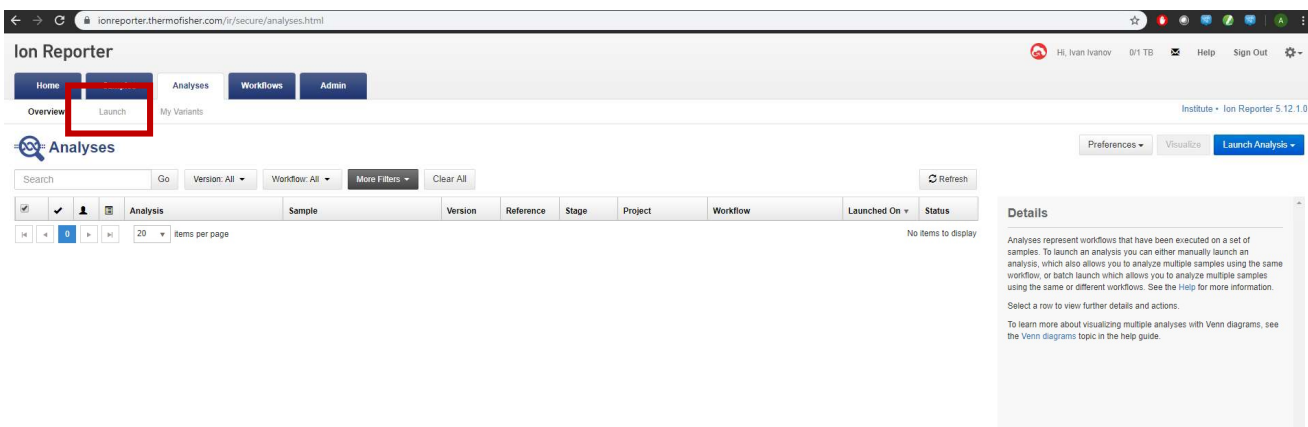

3. В поиске наберите idp и выберите “Ampliseq IDP single sample”. Нажмите “Next”

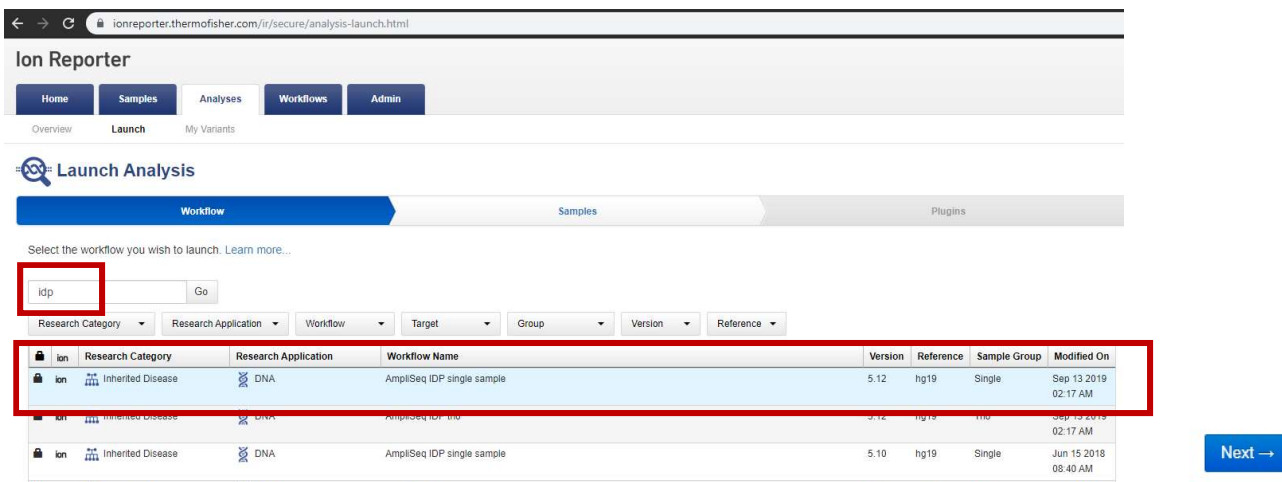

4. В поиске наберите idp и выберите “Demo AmpliSeq IDP Daughter”. Нажмите “Next”

Ion Reporter

Home Samples Analyses Workflows Admin

Overview Launch My Variants

Launch Analysis

Workflow Samples Plugins

Select the sample you wish to analyze. You can select multiple samples and each one will be treated as a separate analysis. [Learn more...](#)

idp Go

Samples

	Sample	Gender	Sample Type	Index	Library ID	Sample ID
<input checked="" type="checkbox"/>	Demo AmpliSeq IDP Daughter	Female	DNA	Unknown	User, Ion	Dec 08 2013 07:39 PM
<input checked="" type="checkbox"/>	Demo AmpliSeq IDP Father	Male	DNA	Unknown	User, Ion	Dec 08 2013 07:39 PM
<input type="checkbox"/>	Demo AmpliSeq IDP Mother	Female	DNA	Unknown	User, Ion	Dec 08 2013 07:39 PM

1 - 3 of 3 items

Previous Cancel Next

5. Нажмите “Next”

Ion Reporter

Home Samples Analyses Workflows Admin

Overview Launch My Variants

Launch Analysis

Workflow Samples Plugins

Select the plugins you wish to include in your analysis. Pricing information for your analysis will be shown on the next page. [Learn more...](#)

In-Analysis Plugins

There are no plugins currently installed. Please contact your organization administrator.

Previous Cancel Next

6. Нажмите “Launch Analysis”

Ion Reporter

Home Samples Analyses Workflows Admin

Overview Launch My Variants

Launch Analysis

Workflow Samples Plugins

Analysis ready to launch!

Review the selected options, name your analysis and then launch it.

Analysis Name: Demo AmpliSeq IDP Mother\_c1529\_2019-10-15-16-41-17-360  
(Demo AmpliSeq IDP Mother)

Description: Optional

Launch Analysis

Previous Cancel

Готовый анализ будет находиться во вкладке Analyses -> Overview